



World's most affordable long duration Pocket Nebulizer

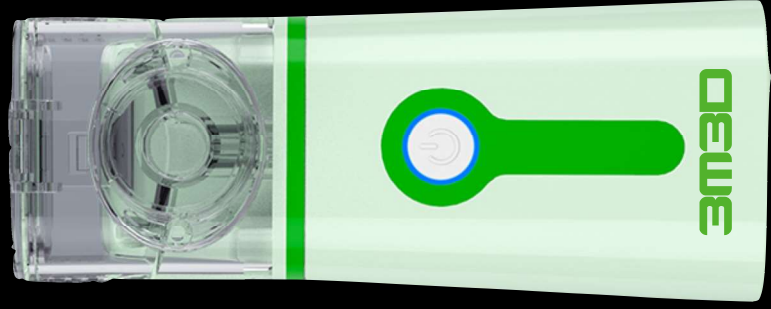

3M3D Pocket NanoNeb is a portable nebulizer that works on 4th generation technology. The nebulizers are able to generate finer nebulized particles and leave low residues as well. The M30 comes with automatic shut of and low battery indication features.

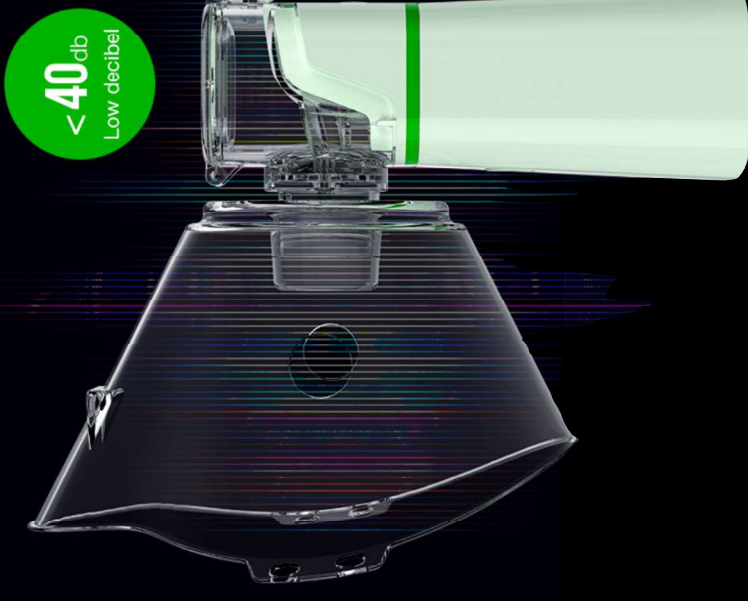

The Particles are exquisite and meticulous.

~3.5 μ m

Median particle size

>1000

Average number of holes in distribution



One-button maintenance in cleaning mode

Self-cleaning function: press for three seconds to start the self-cleaning procedure and maintain without worry.



Automatic shutdown to prevent dry burning

The machine shuts off automatically when lack of water to prevent dry burning and does not need a probe to detect the liquid level.





Automatic shutdown to prevent dry burning

The machine shuts off automatically when lack of water to prevent dry burning and does not need a probe to detect the liquid level.



Small and light to carry freely

Exquisite appearance, small and light, easy to carry.




≤ 0.5 ml
There is little residue
in the solution.


≤ 0.5 ml
There is little residue
in the solution.




Product accessories list

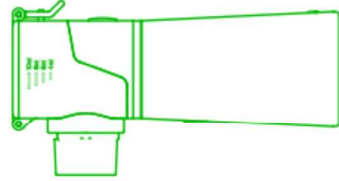

NanoNeb Unit

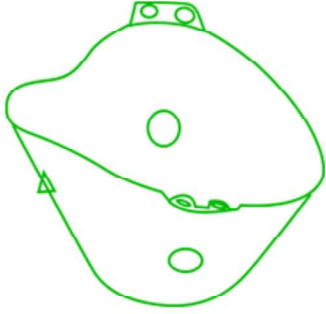


Adult Mask

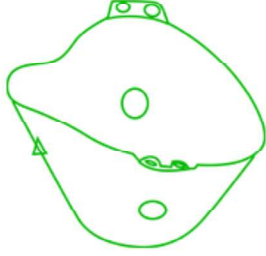

Child Mask



Mouth Piece

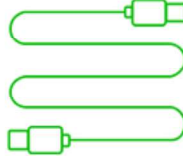

USB Cable



User Manual

3M3D Pocket Nanoneb - Specifications

Device Type	Pocket Nano Nebulizer
Power Supply	Micro USB or 2*AA Batteries
Power Consumption	1.5W
Input	5V/1A
Oscillation frequency	110khz \pm 10%
% of Nebulizer particles with <5μm	100%
Median Particle size	3.5 μ m
Average Nebulization rate	~0.2ml/min
Medication container capacity	Max. 20ml
Weight	115g
Dimensions	45x40x107 (LxWxH)
Operating Temperature range	+5°C ~+45°C
Relative Humidity	<93%RH non-condensing
Storage temperature range	-20°C ~+50°C
Residue	<0.5ml
Noise	<40dB

Model M40

